

GSI 61 USB REMOTE/PRINTER OPTION INSTALLATION INSTRUCTIONS


Items Supplied:


ITEM	QTY	PART #	DESCRIPTION
1	1	512-410500	PCA/GSI61 USBBO ARD
2	1	1761-0224	REMOTE CABLE ASSY
3	1	017-477000	LABEL,GSI 61USB OPTION
4	2	4314-0030	CABLE TIE, PLASTIC
5	4	7044-0024	SCREW, LOCKING,6-32 X.250 LONG
6	1	7044-0037	SCREW, LOCKING,6-32 X.375 LONG
7	1	4320-0400	CABLE CLAMP .312I.D.
8	1	017-463800	LBL S ET, EMC,MAN,GND, CLEAR(LABEL A)
9	1	085-450700	CBL/USBA/B 2M
10	1	1761-9682	REMOTE/PRINTER SOFTWARE UPGRADE KIT
11	1	1761-0189	ADDEN.MANUAL,REMOTE USB INFO
12	1	1761-0650	CD, USB DRIVERS
13	1	1761-9685	GSI61 REMOTE/PRINTER OPTION INSTALLATION SHEET 3 THROUGH 5.

INSTALLATION INSTRUCTIONS:

1


Disconnect power from the GSI 61 and open the unit housing by removing the three screws from the bottom, along the front edge. Save screws for reinstallation. Place unit flat and raise and lock top case into position.

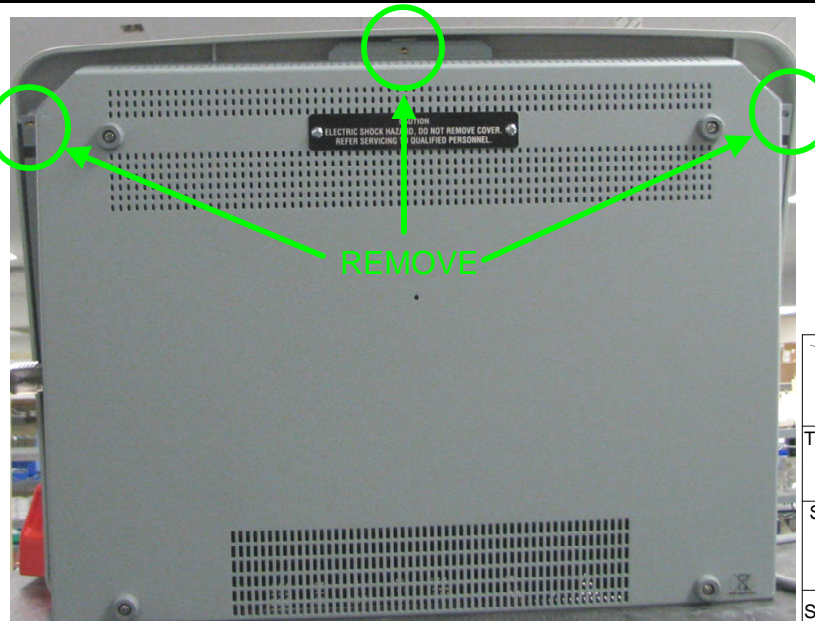

		CARDINAL HEALTH 5225 VERONA RD, BLDG 2 MADISON, WI 53711-4495	
TITLE GSI 61,USB REMOTE/PRNTR OPTION			
SIZE A	DRAWING/PART NO. Reference Only		REV
SCALE/NTS	E-FILE:		SHEET 1 OF 3


2

Remove the blank label on the right side of unit, replace with supplied Remote/Printer Label. **Note:** Remove backing to expose adhesive and use existing screws. Place 017-463800, Label A, above the Remote/Printer label.

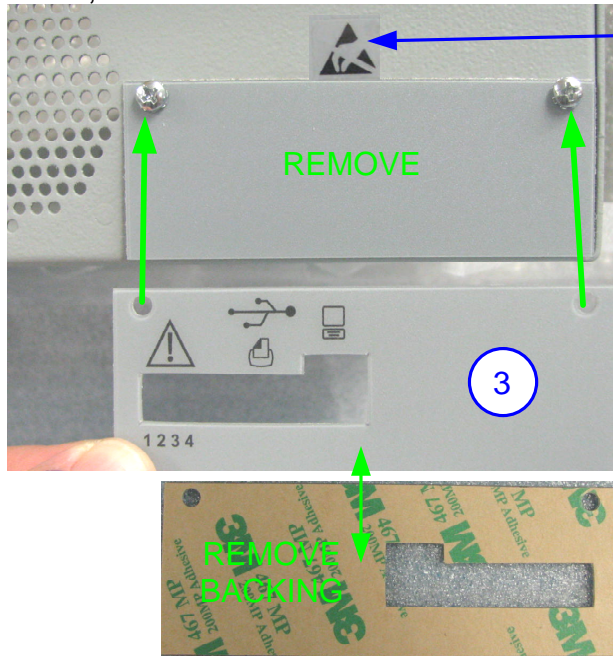


3

Install the Remote/Printer Board using four .250" long screws.

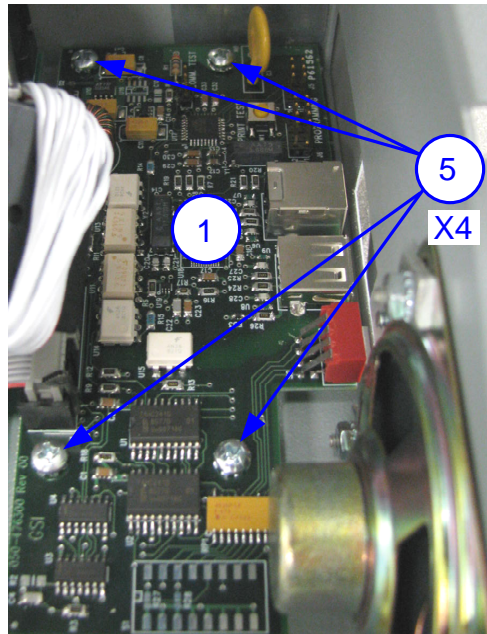


4

Insert and lock Remote Cable Assembly into Remote Board connector J1.

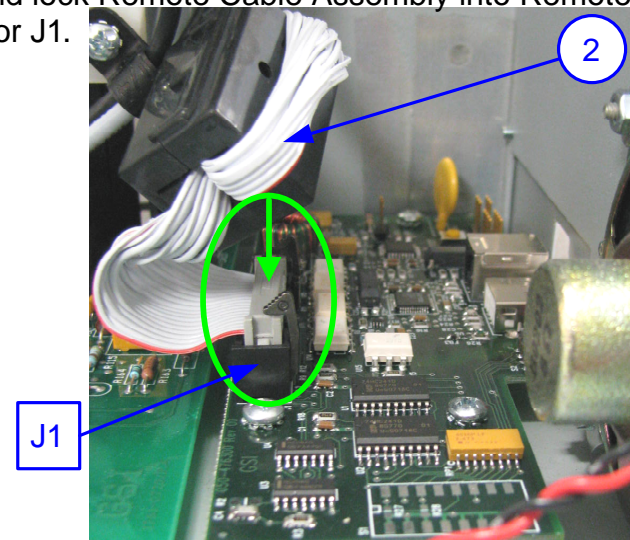

5

Secure Remote Cable to base using cable clamp and one .375" long screw as shown.

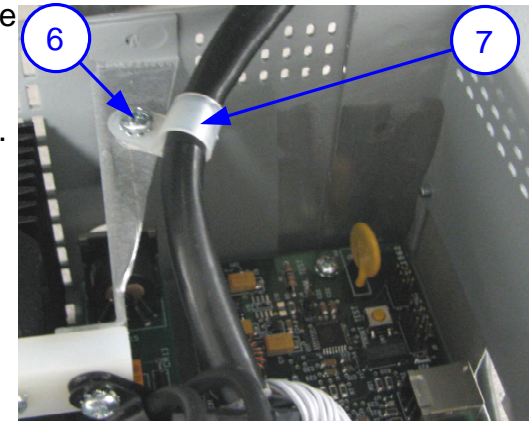

ITEM	QTY	PART #	DESCRIPTION
1	1	512-410500	PCA/GSI61 USB BOARD
2	1	1761-0224	REMOTE CABLE ASSY
3	1	017-477000	LABEL,GSI 61 USB OPTION
5	4	7044-0024	SCREW, LOCKING, 6-32 X .250 LONG
6	1	7044-0037	SCREW, LOCKING, 6-32 X .375 LONG
7	1	4320-0400	CABLE CLAMP .312 I.D.
8	1	017-463800	LBL SET, EMC, MAN, GND, CLEAR (LABELA)
9	1	085-450700	CBL/USB A/B 2M


TITLE
GSI 61,USB REMOTE/PRNTR OPTION

SIZE
A
Reference Only

SCALE:NTS E-FILE: SHEET 2 OF 3

6


Route Remote Cable inside the top case and insert other end of Remote Cable to Main Board, connector J23 and secure cable to top case using two cable ties as shown.

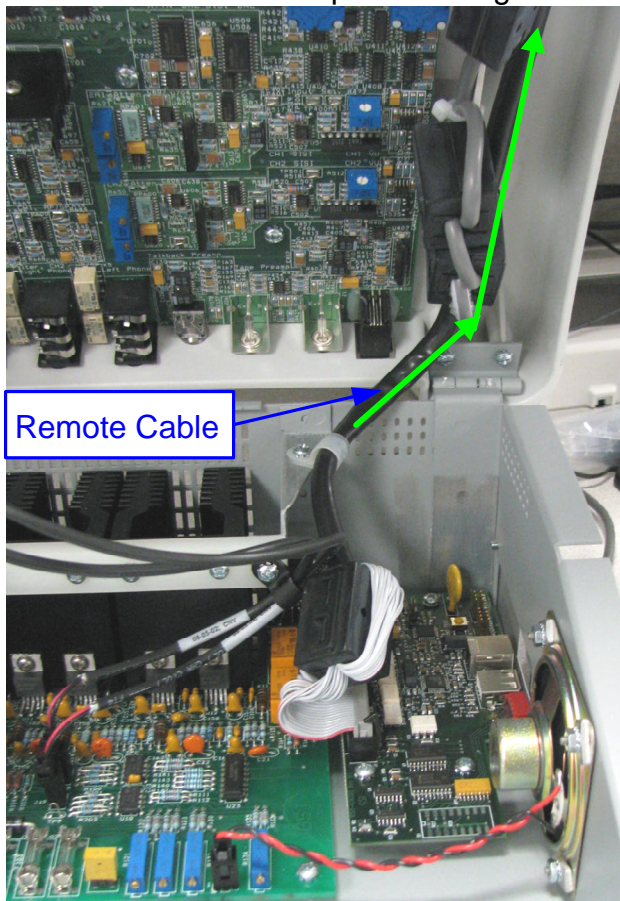

Remote Cable

INSERT AND
LOCK

J23

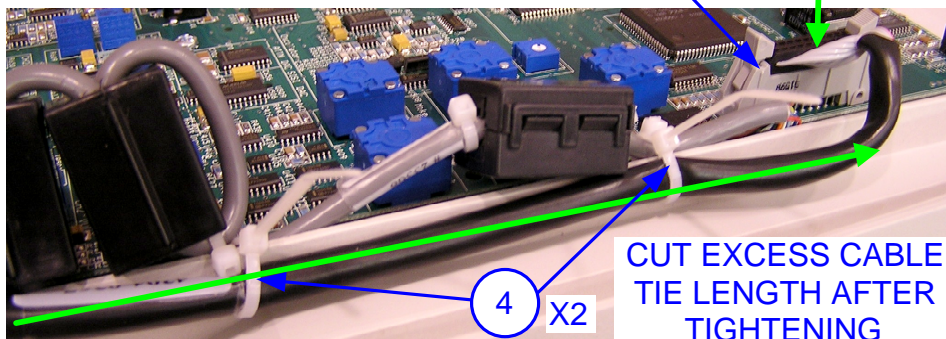

CUT EXCESS CABLE
TIE LENGTH AFTER
TIGHTENING

4 X2

7

Follow instructions on 1761-9682, Remote/Printer Software Upgrade Kit, if necessary for your unit. 10

8

Close the unit housing and replace the three screws removed in step one.

9

Print an Audiogram by pressing, **"DISPLAY AUDIOGRAM"**, **"SAVE"** and **"PRINT AUDIOGRAM"**.

10

See 1761-0189, Addendum, for installing drivers using the supplied 1761-0650, USB Drivers CD, required to connect to remote PC. 11 12

ITEM	QTY	PART #	DESCRIPTION
4	2	4314-0030	CABLE TIE, PLASTIC
10	1	1761-9682	REMOTE/PRINTER SOFTWARE UPGRADEKIT
11	1	1761-0189	ADDEN. MANJAL,REMOTE USB INFO
12	1	1761-0650	CD, USB DRIVERS

TITLE

GSI 61,USB REMOTE/PRNTR OPTION

SIZE

A

Reference Only

REV

SCALE/NTS

E-FILE:

SHEET 3 OF 3